



# WMD

# SWAG



---

Copyright © 2002-2005 by Arkansoft  
All rights reserved.

Sales: [sales@arkansoft.com](mailto:sales@arkansoft.com)  
Technical Support: [support@arkansoft.com](mailto:support@arkansoft.com)

---

## USING WMD SWAG

You can probably use WMD SWAG quite well without reading any of the following information, but if you will take the time to read on, you will discover some helpful information that might make the program easier and more useful to you.

Context-sensitive on-line help is also available in most parts of the program. Just select for “Current screen” on the “About” menu at the top of the screen. (On the handheld computer you need to tap the drop-down menu button to expose the About menu.)



---

## REGISTRATION

The first time you start WMD SWAG for the handheld computer (by tapping its “gas mask” icon on the PDA main menu) you will see a Trial Version screen. When you purchase WMD SWAG you will receive a registration code to convert the trial version to a registered version.

To register the software, tap the [Register] button on the startup screen then use the on-screen number and letter buttons to enter your registration code and submit your registration.

If you enter the code correctly you will get a “registration successful message”.

**TO PURCHASE A REGISTRATION NUMBER**, visit the Arkansoft website at <http://www.arkansoft.com> and click on WMD SWAG in the list of Arkansoft software products.

Then scroll down on the WMD SWAG page and click a logo to get it from PalmGear or Handango. This will take you to a sales and download site.

---

## THE “MAIN SCREEN”

Tap the [Create a New Analysis] button to enter conditions you observe and obtain an educated GUESS of what agent(s) might cause these conditions.



When a new analysis has been created, additional buttons appear on the Main Screen. A [Back] button in the upper right corner of the screen lets you go back to the analysis screens for review or to make changes. A [Next] button lets you go straight to the final analysis screen.



When a new record has not yet been saved, a [Save This Analysis] button appears on the Main Screen to let you, optionally, save the current analysis for future reference.

If you have previously-saved records in the program database, an [Open a Saved Analysis] button appears on the Main Screen to let you access these records.





Finally, if you have BioChem\* installed on your handheld device, a [BioChem] button will appear in the upper right area of the Main Screen to let you move quickly to the BioChem program. If you have the latest BioChem software, a button will also appear on the BioChem Main Screen that will let you easily return to WMD SWAG.



\*BioChem 1stRespondER is a handheld computer program that lists valuable information for paramedics and other disaster first responders to use for WMD data such as first aid advice, recommended drugs, dosages, and lab tests. It's available by visiting <http://www.arkansoft.com>


---

## USING THE DATA ENTRY SCREENS

Just tap the words or checkboxes for any conditions you observe. You can tap again to 'uncheck' a box if you need to change it.

When finished with one screen, tap the  button in the upper right corner of the screen to move to the next screen. Continue in like manner to visit and enter data on each screen available. To go back to previous screens, use the  button.


After the last data entry screen is completed, tapping  will take you to the results screen. You can also, at any point during data entry, go directly to the results screen by tapping the  button.

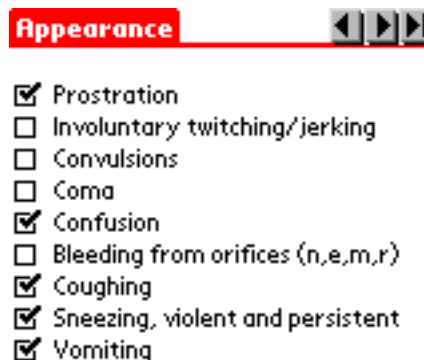


- Prostration
- Involuntary twitching/jerking
- Convulsions
- Coma
- Confusion
- Bleeding from orifices (n,e,m,r)
- Coughing
- Sneezing, violent and persistent
- Vomiting


## USING THE RESULT SCREEN

This screen shows you the program's best, scientific, educated GUESS of what agent or agents are likely suspects to cause the conditions you have entered. The Confidence Level column indicates the degree of likelihood that a particular agent is causing the conditions. The confidence levels may actually be higher than they appear, since several indicators in the environmental section will rarely be seen at all, and most are general indicators that a non-specific WMD event (or industrial accident) has occurred in the vicinity.

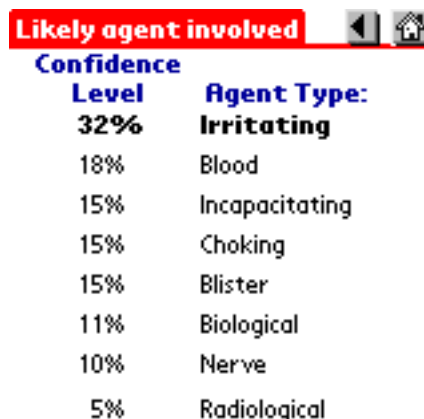


Confidence Level	Agent Type:
32%	Irritating
18%	Blood
15%	Incapacitating
15%	Choking
15%	Blister
11%	Biological
10%	Nerve
5%	Radiological

The  button will let you return to previous screens to observe or change the data entered.

The  button will take you back to the Main Screen.

## SAVING RECORDS

The program lets you save data on an unlimited number of incidents, and to retrieve the incident record. To save a record, first create a new record, checking off any conditions you observe. Then return to the Main Screen and tap the [Save This Analysis] button. You're asked to enter a name for the incident. On the handheld computer, the current system date is filled in for you by default, but you

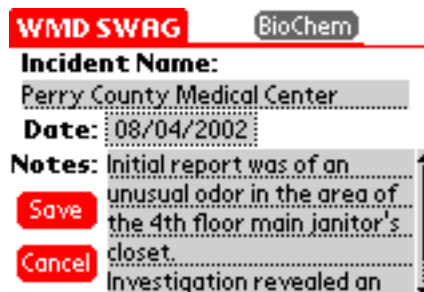

**WMD SWAG** BioChem

**Incident Name:**  
Perry County Medical Center

**Date:** 08/04/2002

**Notes:** Initial report was of an unusual odor in the area of the 4th floor main janitor's closet. Investigation revealed an

Save Cancel

can tap the date and select a different date, if desired. You can also enter notes that will become a permanent part of the record.


---

## THE SAVED RECORDS SCREEN

Once a record has been saved, tapping the [Open a Saved Analysis] button takes you to the Saved Records screen. This screen shows any records that you've previously created and saved in the program database. If more than six records are in the database, scroll arrows appear on the right side of the record list to let you access records that aren't currently displayed.

To select a record on the list, tap its name once. The record's date and notes are displayed. A trashcan button appears. Tapping the trashcan lets you delete the selected record, if desired.

To open a record and review the conditions that were checked when the record was created, tap/click the selected record a second time.

To return to the Main Screen without opening a record, select the  button in the upper right area of the screen.

---

## THE TOP MENU


Some "utility" functions are located on the dropdown "top menu" at the top of the program screen. On the handheld computer, tap the "title bar" at the top of the screen, or tap the PDA's dropdown menu button in the lower left corner of the Graffiti area to open the top menu.

### ***On the "Save Record" Edit menu...***

Copy and Paste can allow you to copy information in another WMD SWAG record, or anywhere else on the PDA, to the memory "clipboard" then paste it into a new WMD SWAG record. Or you can copy information in a record and paste it into a PDA e-mail, the PDA memo pad, or anywhere else you may need it.

To copy text, drag the PDA's stylus over the text to "highlight" it, then dropdown the edit menu and select copy.

To paste text, tap the stylus where you want to insert the copied text then select Paste from the dropdown menu.

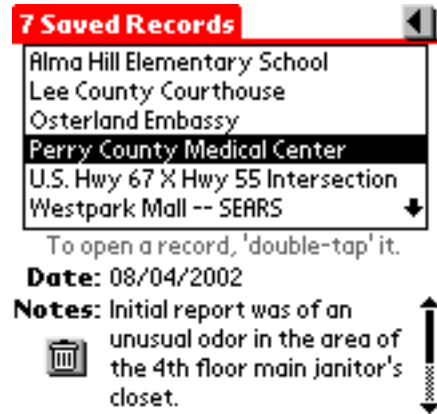

### ***On the handheld Beam menu (Windows Print menu)...***

The handheld Beam menu lets you beam a trial version of WMD SWAG to another handheld computer. You can also, optionally, beam your handheld's incident data files to another device.

On incident record screens, you can generate a printed report of the current record. On the handheld computer, the print command is on the Beam menu. If you have an IR-capable printer or serial cable connection and PalmPrint or SCSPrint on your device, you can beam reports directly to a printer. A free trial copy of SCSPrint is included with WMD SWAG so you can give it a try.

In windows, the print option will export the report as a DOC text file suitable for opening/editing/printing as needed from any Windows word processor. You are given the opportunity to name the DOC file and specify its storage location on your hard disk, so you can easily find it later to print it.

If you try to print when you don't currently have a record's data loaded or entered into memory, you will get a "Nothing to print" message. To print, you must first create, name, and save a new record, or open a previously saved record.

There's no beaming available in the Windows version, of course.

### ***On the About menu...***

You can get information on how to use the currently-viewed screen. You can also get the version number of WMD SWAG or get information on how to contact Arkansoft.

## REMOVING WMD SWAG FROM YOUR DEVICE

Should you need to remove WMD SWAG from a handheld computer, just use the standard Palm OS Delete function. Tap the drop-down menu button from the device “main menu” and select Delete. Scroll down and highlight WMD SWAG on your program list and tap the Delete button.

To remove the WMD SWAG Windows version, first go to the WMD SWAG program group in the Windows Start>Programs section and select “d) Remove WMD SWAG (ARKW) Conduit”, then confirm that you want to remove the program conduit. Next, return to the WMD SWAG Windows version, first go to the WMD SWAG program group in the Windows Start>Programs section and select Uninstall to remove the Windows version of WMD SWAG.